



PAL-BX/RI features a dual scale. Brix (wide range 0.0-93.0%) and Refractive Index (resolution 0.0001).

Digital Hand-held "Pocket" Dual Scale Refractometer

PAL-BX/RI

Cat.No.3851

Measurement Method



1. Place the sample onto the prism surface.



2. Press the START key.




3. The measurement value will be displayed in 3 seconds.

4. Press the START key for 2 seconds while the Brix (%) is indicated, to display the Refractive Index and temperature. Press the START key (one time) to return to the Brix (%) scale.

Example of display

Display of the measurement value "1.3425" and "20°C".

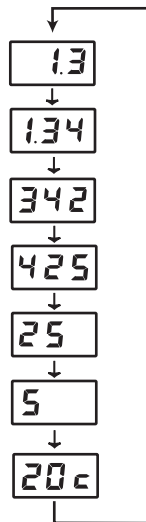

The measurement value for Refractive Index will scroll from right to left.



Quick & Easy cleanup!



Actual size

Features

- Requires only a few drops of sample
- Measurement in 3 seconds
- Ergonomically designed for easy, one-handed operation
- Light & Compact, 100g
- Designed with a strap holder option
- Highly water resistant (IP65)
- Convenient storage case
- Battery power indicator

Specifications

Measurement range	Brix : 0.0~93.0%	Refractive Index (RI) : 1.3306~1.5284
Resolution	Brix : 0.1%	Refractive Index : 0.0001
Measurement accuracy	Brix : ±0.1%	Refractive Index : ±0.0003 (Water at 20°C)
Measurement temperature (Automatic Temperature Compensation)	10~75°C (Brix only)	
Ambient temperature	10~40°C	
Power supply	AAA alkaline battery × 2	
Battery life	Approx. 11, 000 measurements (using new alkaline batteries)	
International Protection Class	IP65 Water resistant	
Dimensions and weight	55(W) × 31(D) × 109(H)mm, 100g (main unit only)	

* The displayed measurement value for Refractive Index is not temperature compensated. The reading is the Refractive Index of the sample at the temperature displayed.

All ATAGO refractometers are designed and manufactured in Japan.



ATAGO products are acknowledged by HACCP, G.M.P. and G.L.P.

ATAGO CO.,LTD.

Headquarters: 32-10, Honcho, Itabashi-ku, Tokyo 173-0001 Japan
 TEL: 81-3-3964-6156 FAX: 81-3-3964-6137
 overseas@atago.net http://www.atago.net/

ATAGO U.S.A.,Inc.

12011 NE First Street, Bldg. C, Suite 110 Bellevue, WA 98005 U.S.A.
 TEL: 1-425-637-2107 FAX: 1-425-637-2110
 customerservice@atago-usa.com http://www.atago.net/

ATAGO INDIA Instruments Pvt. Ltd.

A-501, Mangal Aarambh Business Centre, Near Kora Kendra, Off. S.V. Road, Borivali (West), Mumbai-400 092 India
 TEL: 91-22-2833-8038 / 8076 FAX: 91-22-2899-8164
 customerservice@atago-india.com



※Specifications and appearance are subject to change without notice.